




LAYERING TECHNIQUE

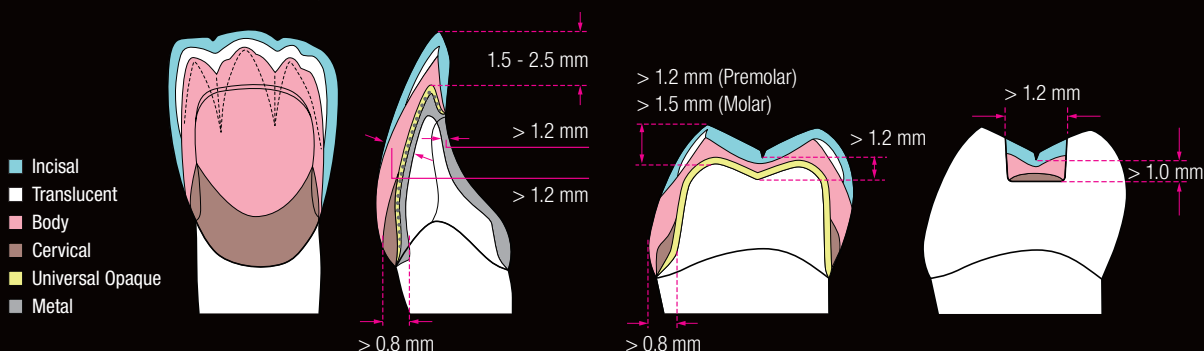

ENAMEL COLOUR IMAGE

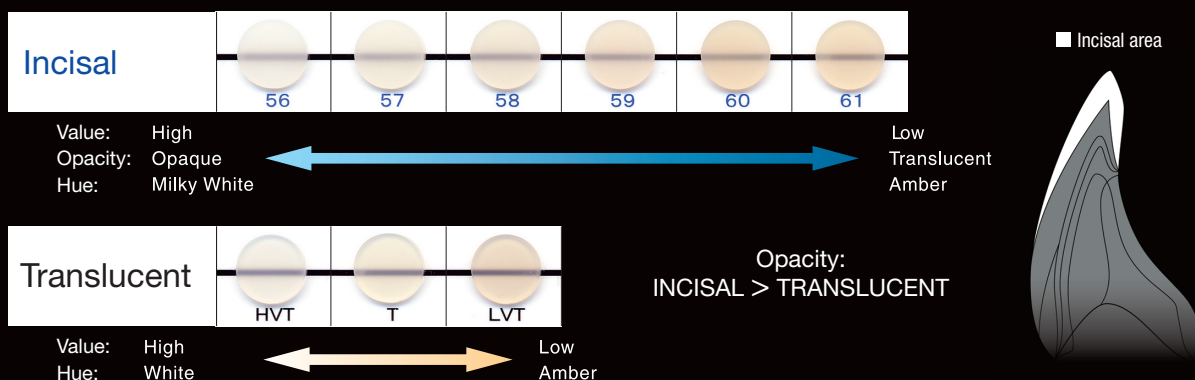

LAYERING TECHNIQUE FOR GINGIVAL REPRODUCTION

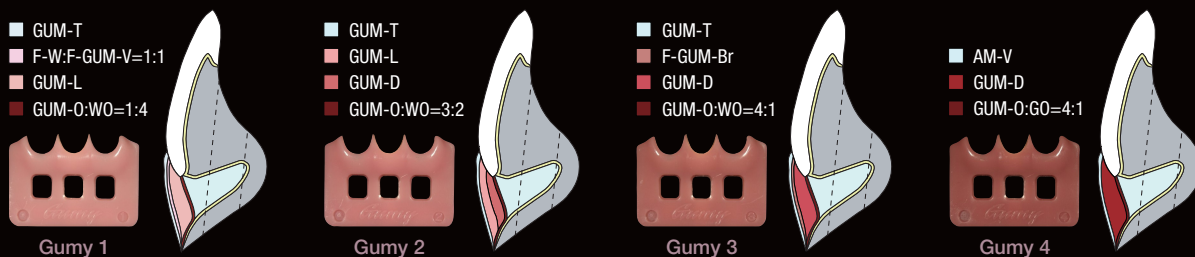

LIGHT-CURING TIME

		Solidilte LED	Solidilte V
SHOFU Universal Opaque	Pre-Opaque	30	60
	Opaque	90	180
Prepolymerisation		30	60
Final polymerisation		90	180

(seconds)

COLOUR	W0	W1	W2	W3	A1	A2	A3	A3.5	A4	rootA
UNIVERSAL OPAQUE	W0	W1	W2	W3	A1	A2	A3	A3.5	A4	rootA
CERVICAL						AC1		AC2		
OPAQUE DENTIN					ODA1	ODA2	ODA3	ODA3.5	ODA4	ODrootA
BODY	W0B	W1B	W2B	W3B	A1B	A2B	A3B	A3.5B	A4B	rootAB
INCISAL		56		57		58		59		60

COLOUR	B1	B2	B3	B4	C1	C2	C3	C4	D2	D3	D4
UNIVERSAL OPAQUE	B1	B2	B3	B4	C1	C2	C3	C4	D2	D3	D4
CERVICAL		BC1		BC2		CC1		CC2		DC1	DC2
OPAQUE DENTIN	ODB1	ODB2	ODB3	ODB4	ODC1	ODC2	ODC3	ODC4	ODD2	ODD3	ODD4
BODY	B1B	B2B	B3B	B4B	C1B	C2B	C3B	C4B	D2B	D3B	D4B
INCISAL	57	58	59	60	58	59	60	59	60	59	59

The printed colours might be different from the actual colours in this system.

COLOUR	R2	R3	R3.5
UNIVERSAL OPAQUE	R2	R3	R3.5
OPAQUE DENTIN	ODR2	ODR3	ODR3.5
BODY	R2B	R3B	R3.5B
INCISAL	58	59	

UNIVERSAL OPAQUE											
TRANSLUCENT											
INCISAL											
CONCENTRATE											
FLOWABLE COMPOSITE RESIN											
											